

<b>Model Number</b>	BA401
<b>Weight Capacity</b>	450 lbs.
<b>Maximum Speed<sup>1</sup></b>	Up to 11 mph at 450 lbs. Up to 11.4 mph at 200 lbs.
<b>Ground Clearance<sup>3</sup></b>	8.25" at mid-frame 5.5" at motor
<b>Turning Radius<sup>3</sup></b>	94.75"
<b>Overall Length<sup>3</sup></b>	61.25"
<b>Overall Width<sup>3</sup></b>	32.5"
<b>Seat-to-Ground Height<sup>3</sup></b>	29-30.5"
<b>Seat-to-Deck Height<sup>3</sup></b>	20.5-22"
<b>Front Tires<sup>3</sup></b>	14.5" pneumatic
<b>Rear Tires<sup>3</sup></b>	14.5" pneumatic
<b>Suspension</b>	Front and rear (CTS Suspension)
<b>Range Per Charge<sup>1,2</sup></b>	75 AH: Up to 13.2 miles at 450 lbs. Up to 20.2 miles at 200 lbs. 100 AH: Up to 17.4 miles at 450 lbs. Up to 26.6 miles at 200 lbs
<b>Total weight<sup>5</sup></b>	274 lbs.
<b>Standard Seating<sup>3</sup></b>	Black vinyl high-back with limited recline memory foam seat with sliders Weight: 48.6 lbs. Size: 16"x16-18" or 18"x18-20"
<b>Drive System</b>	24-volt DC motor, rear-wheel drive, dual in-line
<b>Braking System</b>	Regenerative & electromechanical
<b>Battery Requirements<sup>4,6</sup></b>	(2) 75 AH (100 AH optional) Weight: 48 lbs. (75 AH) 60.5 lbs. (100 AH)
<b>Battery Charger</b>	8-amp, off-board
<b>Warranty</b>	Frame: lifetime limited Drivetrain: 2-year limited Electronics: 2-year limited Battery: 6-months

**\*HCPCS Code: K0899**

The Baja® Wrangler® 2 is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

- 1 Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition, and tire condition. This specification can be subject to a variance of (+ or -) 10%.
- 2 Tested in accordance with ANSI/RESNA, MC Vol2, Section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- 3 Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or -) 3%.
- 4 Battery weight may vary based on battery manufacturer.
- 5 Includes standard seat and batteries.
- 6 AGM or gel cell type required.



# BAJA® WRANGLER® 2

SHOWN WITH OPTIONAL  
CANOPY UPGRADE



800-800-8586 • [bajamobility.com](http://bajamobility.com)



HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

FDA Class II  
Medical Device

## THE SPIRIT OF ADVENTURE

WRANGLER 2 IS MADE IN  
**AMERICA**  
OF FOREIGN AND DOMESTIC PARTS



# BAJA® WRANGLER® 2

OPTIONAL hard top canopy with durable, powder-coated steel frame, polycarbonate windshield, panoramic rearview mirror, and rear mounting brackets

Limited recline memory foam seat with manual adjustment, sliders, and height, width, depth, and angle adjustable armrests

Heavy-duty seat post

Rear backup sensor with LED lights

FDA Class II Medical Device\*



Rearview mirrors (standard)

User friendly console with tiller adjustment lever

Hydraulic hand brake

USB charger in tiller

Storage compartment

Front LED headlights, directional signals and hazard lights

Front and rear CTS suspension and additional steering dampening


## REINVENTING THE SPIRIT OF ADVENTURE

Rugged and built to last, the Baja® Wrangler® 2 is everything you want in an outdoor scooter. Equipped with front and rear suspension, this aggressive scooter is designed for cruising on walking trails or driving around your property. Make your adventures legendary with the Baja Wrangler 2!

## FEATURES

- Industry-leading front and rear suspension
- Steering dampening system absorbs rugged, uncontrolled movement, reducing energy required to control handles
- Choice of 16"x16-18" or 18"x18-20" black vinyl, limited recline memory foam seat with manual adjustment, heavy-duty seat post, sliders and height, width, depth, and angle adjustable armrests (optional seat sizes of 20"x18-20" and 22"x20-22" available)
- User friendly LED console with variable speed control displays time, temperature, and mileage, and also includes horn and light buttons
- Adjustable delta tiller with ergonomic, wraparound handles, XLR charging port and USB mobile device charger
- Easily accessible tie-down points (for safe transport of unoccupied scooter)
- LED headlights, hazard lights, and directional signals
- Rear backup sensor LED lights provide safer reversing

## COLORS OPTIONS

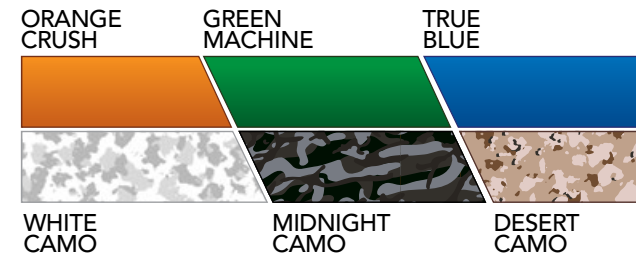

## OPTIONAL ACCESSORIES

- 3" single or double cane/crutch holder
- 50/60/70" lap belt
- Cup holder
- Cell phone holder
- Front basket
- Oxygen holder
- Rear basket
- Saddlebag
- Storage pod
- Walker holder

### NEW Hard top canopy and optional canopy accessories:


- Tool clips (pair)
- Jumbo rear basket
- Canopy weather cover

### NEW Depth-adjustable high-back seat:

- 32"W x 18-20"D (without sliders)



User-friendly console

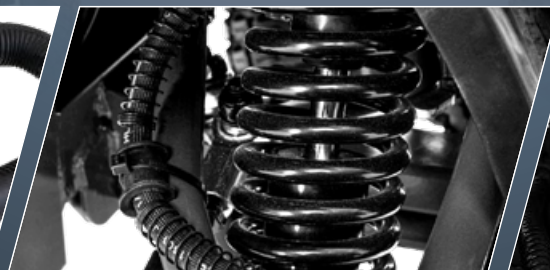

State-of-the-art suspension



Full LED lighting package