



"THIS WHEELCHAIR HAS GIVEN ME BACK MY FREEDOM AND OPENED UP MY WHOLE LIFE AGAIN."

Lisa
New Zealand

MAGIC MOBILITY FRONTIER V6 ALL TERRAIN

**AUSTRALIAN
MADE**
-
BUILT TO LAST
-
**CUSTOMISED
TO YOUR
REQUIREMENTS**



FULL LIGHTING PACKAGES

Been out having too much fun and it's getting late? Put your lights on. Fitted on rubber brackets.



UP CLOSE AND PERSONAL

Side transfers aren't the only reason you need to flip your armrest out of the way.



REHAB SEATING

Infinite options for infinite body shapes. Our rehab seating range is easily adjustable or can be custom designed to a fixed position. It's whatever you need.

TAKE CONTROL

A full range of advanced control systems and environmental controls are available. Easy to program and easy to use.



MAGIC MOBILITY FRONTIER V6 ALL TERRAIN



Look no further! This is it, our most popular wheelchair worldwide due to its versatility. An engineering beauty that adapts to any environment. The Frontier V6 All-Terrain is the ideal outdoor wheelchair, yet its size makes it a capable indoor performer. Its suspension and outdoor performance lets you tackle the bumpy terrain life throws at you.

Basic specifications:

- Length: 1050mm
- Width: 710mm
- Seat to floor height: 435mm
- Turning radius: 550mm
- Weight limit: 182kg (400lbs), 155kg (340lbs) with lift or tilt
- Top speed: 10kph (6.2mph)
- Motors: 800W Gear in Line
- Wheels: 14" (145/70-6) black knobby
- Castors: 210x60mm (2.50-4) black
- Batteries: 73Ahr Gel
- Penny and Giles controller with LCD display (R-Net)
- Centre post legrest with one piece footplate
- Sling backrest and rehab seat pan (from 12"-23" wide)*
- Crash test compliant tie down points.

Rehab positioning:

- Power options: 300mm seat lift, tilt (45 degrees), power recline, lie down, power elevating legrests, transfer tilt (8 degrees and up to 70mm floor to seat height reduction)
- Seating and Backrests: MPS, rehab backrest and adjustable seat pans, foam in place, custom seating (built by our in house team of fabricators), a choice of off the shelf backrests and seat cushions
- Headrests and mounts: MPS, rehab (single and multi axis stem, pivoting, Y-type, L-Type), any custom requirement
- Legrests: Centre post, side mount, swingaway, power elevating
- Footplates: 1 or 2 piece fixed/flip up/angle adjustable, any custom requirement
- A full range of armpads, straps, pads, harness, restraints, supports are all available from our extensive catalogue.

Options and accessories:

- 24V LED Lighting kit
- USB charger
- Jack
- Luggage rack, bag hooks, bags, pockets, keep safe secret boxes, tables, drinks holder, oxy and vent carriers
- Special controls: Attendant, head, chin, foot, sip and puff, environmental controls all available
- Indoor wheel kit: narrower, grey drive wheels and castors, when all-terrain capabilities are not needed. This also reduces the size of the chair further.



COMPACT 40



COMPACT 73



HYBRID

Specifications subject to change. Check online for current version

"WE'VE BEEN
CREATING
CUSTOM MADE
WHEELCHAIRS
FOR 20 YEARS,
SO THAT YOU CAN
LIVE YOUR LIFE"

Ashley Daff
Director




ATTENTION TO DETAIL

Working together, we can figure out how to make your wheelchair uniquely yours. We know the small touches make a big difference.



BUILT TO LAST

We never compromise on quality. With your wheelchair as your mobility, reliability is paramount. We'd much rather spend time improving our products than repairing them.

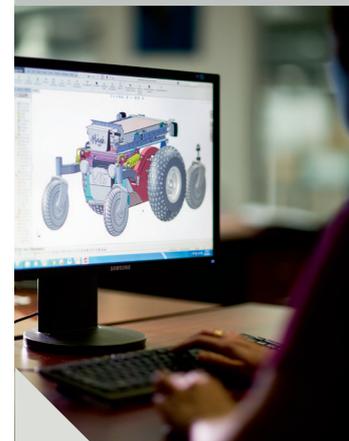

ENGINEERED TO THE LAST BOLT

If we think we can make it better, we will. Our in house engineering team are hands on, working side by side with customers and therapists. The entire Magic team work continuously on how we can improve the experience from choosing a Magic Mobility wheelchair to its performance in the places you want to go.



PROUDLY AUSTRALIAN OWNED, DESIGNED AND MADE

We love that our product is made in Australia. We love meeting our customers, hearing their ever changing needs and taking it back to the drawing board.

*Various model options excluded



Visit us online to see how
we change lives every day:
magicmobility.com.au

Postal Address

Magic Mobility Pty Ltd
2/112 Browns Road
Noble Park VIC 3174 Australia

Tel +61 3 8791 5600

Fax +61 3 9701 0144

Email sales@magicmobility.com.au



Join the
Magic Mobility
Facebook and
LinkedIn
communities!

